



NEW
3rd Edition

MAG-Expert®

The new benchmark in magnetotherapy
in a contemporary design



PHYSIOMED®
technology for therapy

MAG-Expert®

Professional magnetotherapy device for hospitals, rehab and physiotherapy centres



With the **PHYSIOMED App** treatment protocols can be created, archived and processed!



MAG-Expert is available in the following configurations:

MAG-Expert with anthrazit Coil Ø 60 cm integrated in a therapy couch – pillow included – for magnetic field therapy of the spine, pelvis and the whole body



Cylinder, Ø 30 cm for magnetic field therapy of the limbs and head



Magnetotherapy with set of applicators


Two-channel magnetotherapy

MAG-Expert provides magnetic field treatment with a fieldstrength of 1–100 Gauss (adjustable in steps of one Gauss) and a frequency range from 1–100 Hz, with two completely independent channels and a treatment timer. Cylinders of different sizes as well as high performance applicators are available to maximise the efficiency of treatment

MFC technology focuses the magnetic fields almost entirely on the inside of the cylinder. This avoids unnecessary exposure of the treatment team.



The detailed indication menu with practical filter functions

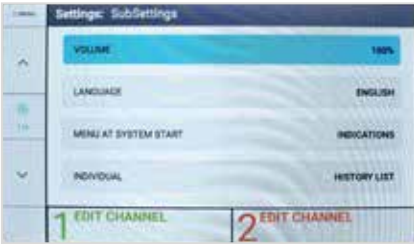

Two-channel magnetotherapy with visible parameters that can be changed at any time

Technical Data	
Protection class	1, Type BF
Power connection	230 V ±10 % or 115 V ±10 %
Mains frequency	50–60 Hz
Fuses	4 A (at 230 V) or 6.3 A (at 115 V)
Power consumption	410 VA
Therapy frequency	1–100 Hz
Power output	1– 00 Gauss
Dimensions (W x H x D)	24,5 x 10 x 31 cm
Weight	6.8 kg
Therapy couch:	
Dimensions (W x H x D)	190 x 75 x 70 cm
Weight	50 kg

Special Features

Magnetotherapy
Magnetotherapy from 1–100 Gauss
Therapy frequency from 1–100 Hz
Two independent channels
MFC technology for focussing magnetic fields
Synchronous mode emission Ch1 and Ch2
Therapy time 0–60 Min
Automatic detection of connected accessories
Therapy couch
Extremely stable treatment couch with pillow
Precisely manufactured sliding for easy and comfortable positioning of the coil

General Features

7" colour touchscreen monitor visualizing all main parameters of the active channels
12 Favorites menu with speed-dial memory for individual device functions
Comprehensive overview of the therapy parameters including all therapy timers
Fastest therapy start: directly, through channel selection, from program memory of indications index, individual or favourites menu
Treatment index with intelligent filtering functions based on body region, sorted in alphabetical order
Extensive memory capacities, using the Individual, history and favourites menu
LED for visual indication of active channel
Audible signal for indication of end of treatment
PHYSIOMED App compatibility
Report generation and error interpretation through QR-Code

Standard accessories

1 Operating instruction
1 Mains cable
1 Test metal ring
1 Connection cable

Accessories depending on the configuration

1 Therapy couch with Coil Ø 60 cm, 230 V/115 V
1 Coil Ø 30 cm, 230 V or 115 V
1 Set of applicators 230 V/115 V for local application



Magnetotherapy

With therapeutic magnetic field devices, pulsating magnetic fields (PEMF) can be applied at various frequencies and intensities to produce good effects, specially to regenerate tissues and bones. Coils and applicators of different sizes can be used.

Magnetic field therapy can produce the following effects:

- Pain relief
- Suppression of inflammation
- Acceleration of healing processes
- Stimulation of blood circulation

Main indications

- **Recent fractures and consolidation delays**
Magnetic fields accelerate the consolidation time in a very high number of cases by stimulating the osteogenic activity at the fracture rime level.

The recovery process of bone discontinuity takes place through the induction of the piezoelectric effect in the connective structures improving the local circulatory conditions.

- **Pseudoarthrosis**
Magnetic fields have been shown to be effective in treating pseudoarthrosis.
- **Osteoporosis**
Magnetic therapy was effective in 70 % of cases with osteoporosis.
- **Inflammatory and degenerative arthropathy**
Both inflammatory and degenerative arthropathies are treated by magnetic therapy. Positive results have been obtained for rheumatoid arthritis in hands and knees, spondylitis ankylopoetica, gonarthrosis and lumbar arthrosis

ADDRESS

PHYSIOMED ELEKTROMEDIZIN AG
Hutweide 10
91220 Schnaittach
Germany

PHONE

+49 9126 2587 - 0

FAX

+49 9126 2587 - 25

E-MAIL

info@physiomed-group.com

WEB

www.physiomed-group.com

